

UNLOADER ONE® X

Unloader One X improves mobility by providing pain relief from knee osteoarthritis and degenerative meniscal tears. It features Össur's unique patented and clinically proven 3-point Leverage System that generates unicompartmental unloading of the knee joint.



OA KNEE


INDICATIONS FOR USE

- Moderate to severe unicompartmental knee osteoarthritis
- Degenerative meniscal tears
- Other knee conditions that may benefit from unicompartmental unloading

FEATURES AND BENEFITS

- 3-Point Leverage System consisting of two Dynamic Force System (DFS) Straps, rigid shells and a single upright hinge contribute to a unicompartmental unloading of the knee joint.
- The unloading is effective only when needed, no condyle pressure during sitting or swing phase.
- A slender perforated frame with soft edges, contours easily to individual anatomies and remains discreet whether worn over or under clothing.
- Skin friendly liners and pads as well as a smooth inner surface of shells reduce probability of pressure points.
- The gastroc strap is positioned over thickest part of calf to help prevent migration.
- Tool-free SmartDosing® dials enable the patient to be in control of their pain level by fine-tuning the strap tension.
- Color coded touchpoints and Quick-Fit buckles increase user-friendliness.
- Alignment guide on hinge facilitates correct positioning on the leg.
- Entire brace is machine washable (washbag included with product).
- Easily adaptable to individual needs and anatomical shapes.
- Field serviceable.

REQUIRED MEASUREMENTS

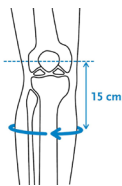

Measure calf circumference
15 cm (6") below mid-patella.
Leg in extension.

Size	Calf circumference 15 cm (6") below mid patella	
	cm	inch
XS	25 - 31	9¾ - 12¼
S	31 - 36	12¼ - 14
M	36 - 43	14 - 16¾
L	43 - 51	16¾ - 20
XL	51 - 62	20 - 24¾

UNLOADER ONE® X

ORDERING INFORMATION

STANDARD

Size	Right Medial	Left Medial	Right Lateral	Left Lateral
XS	UOX0518711	UOX0618711	UOX0528711	UOX0628711
S	UOX0518712	UOX0618712	UOX0528712	UOX0628712
M	UOX0518713	UOX0618713	UOX0528713	UOX0628713
L	UOX0518714	UOX0618714	UOX0528714	UOX0628714
XL	UOX0518715	UOX0618715	UOX0528715	UOX0628715

Recommended for patients between 1.65 and 1.92 meters tall (approx. 5'5" - 6'4"). Hinge is placed on the side of the affected compartment.

SHORT

Size	Right Medial	Left Medial	Right Lateral	Left Lateral
XS	UOX1518711	UOX1618711	UOX1528711	UOX1628711
S	UOX1518712	UOX1618712	UOX1528712	UOX1628712
M	UOX1518713	UOX1618713	UOX1528713	UOX1628713
L	UOX1518714	UOX1618714	UOX1528714	UOX1628714
XL	UOX1518715	UOX1618715	UOX1528715	UOX1628715

Recommended for patients under 1.65 meters tall (approx. 5'5"). Hinge is placed on the side of the affected compartment.

ACCESSORIES

EASE-OF-USE KIT

Description	Part#
Kit	UOX800001

Sundial
Buckle finger loop

THIGH EXTENSION LINER, GRAY

Right Medial, Left Lateral	Right Lateral, Left Medial
UOX546101	UOX545101

ACTIVEGRIP CALF LINER, GRAY

Size	Right Medial, Left Lateral	Right Lateral, Left Medial
XS-S	UOX516153	UOX515153
M-XL	UOX516163	UOX515163

DOESKIN CALF LINER, GRAY

Size	Right Medial, Left Lateral	Right Lateral, Left Medial
XS-S	UOX516151	UOX515151
M-XL	UOX516161	UOX515161

SLEEVE

Size	Part#
XS	UOX60001
S	UOX60002
M	UOX60003
L	UOX60004
XL	UOX60005
XXL	UOX60006

UNLOADER ONE® X

Description	Part#
Hinge pad, gray	UOX700011
Extension stop kit	B-301040047
Sensil strap pad, gray (10 pcs)	UOX700193
Sundial	UOX1200001

REPLACEMENT PARTS

Description	Part#
Popliteal pad, gray	UOX700092
Alligator clip for DFS, gray	UOX900002
SmartDosing Dial, gray (1pc)	UOX1000002

UPPER DFS AND CALF STRAP WITH BUCKLE, GRAY

Right Medial, Left Lateral	Right Lateral, Left Medial
UOX406101	UOX405101

Calf & DFS strap with buckle

LOWER DFS AND THIGH STRAP WITH BUCKLE, GRAY

Right Medial, Left Lateral	Right Lateral, Left Medial
UOX406102	UOX405102

Thigh & DFS strap with buckle

SENSIL CALF LINER, GRAY

Size	Right Medial, Left Lateral	Right Lateral, Left Medial
XS-S	UOX516152	UOX515152
M-XL	UOX516162	UOX515162

DOESKIN THIGH LINER, GRAY

Right Medial, Left Lateral	Right Lateral, Left Medial
UOX526101	UOX525101